


Malaffi Provider Portal Quick Reference Guide - Malaffi Image Viewer


The Malaffi Image Viewer displays Radiology Images available for your patients. This ensures that you have a holistic view of your patient record to be able to review and compare images from different visits and different timeframes.

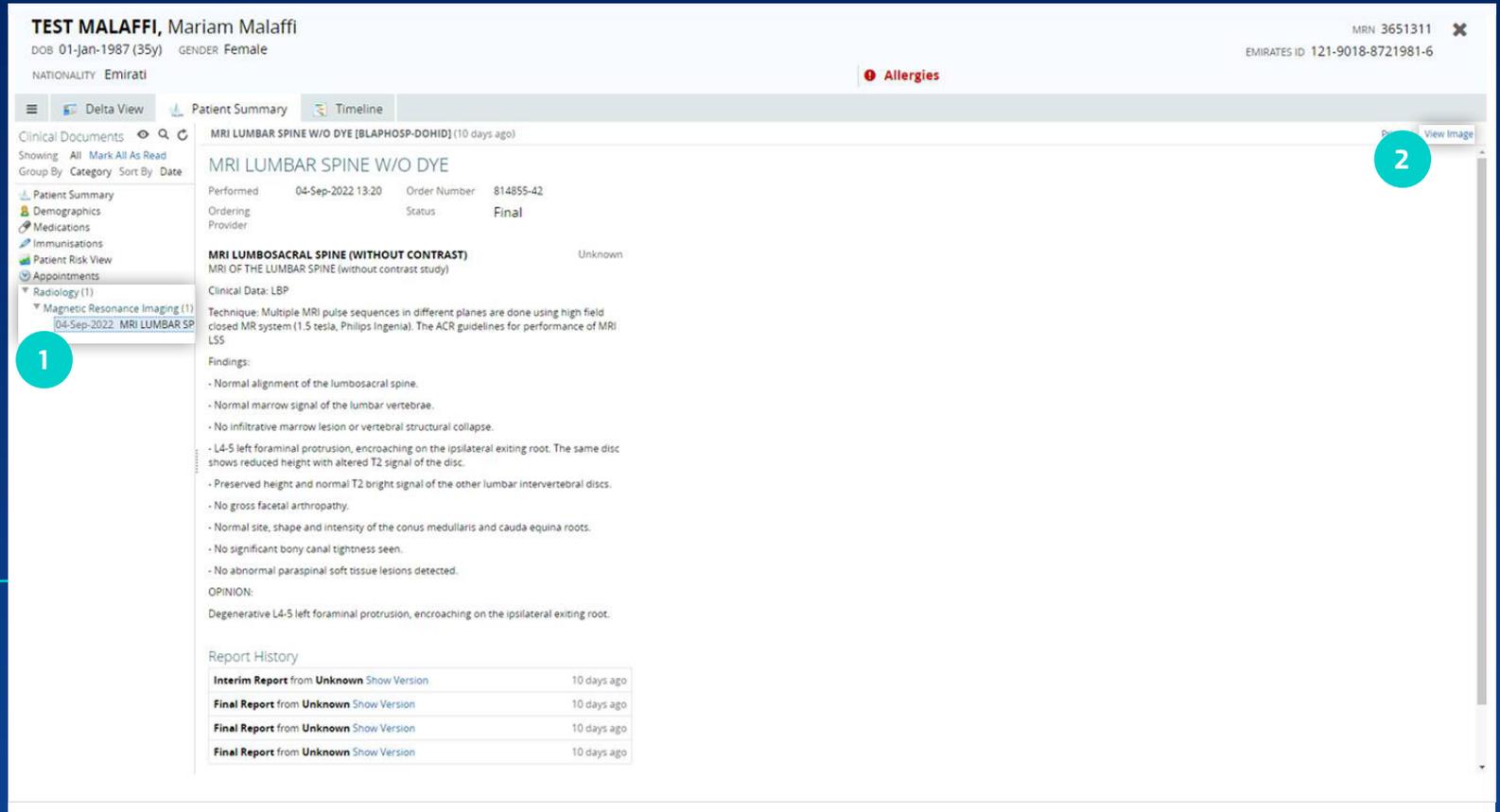

1

Select Radiology link

To access the images, click and expand the *Radiology* folder in the CDV tree. Select the relevant record to view report.

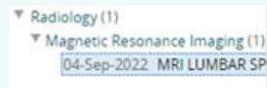


2

'View image' link

If there is an image associated with that report, the *View Image* link will be present. If there is no image associated with the report, the link will not be present.

View Image

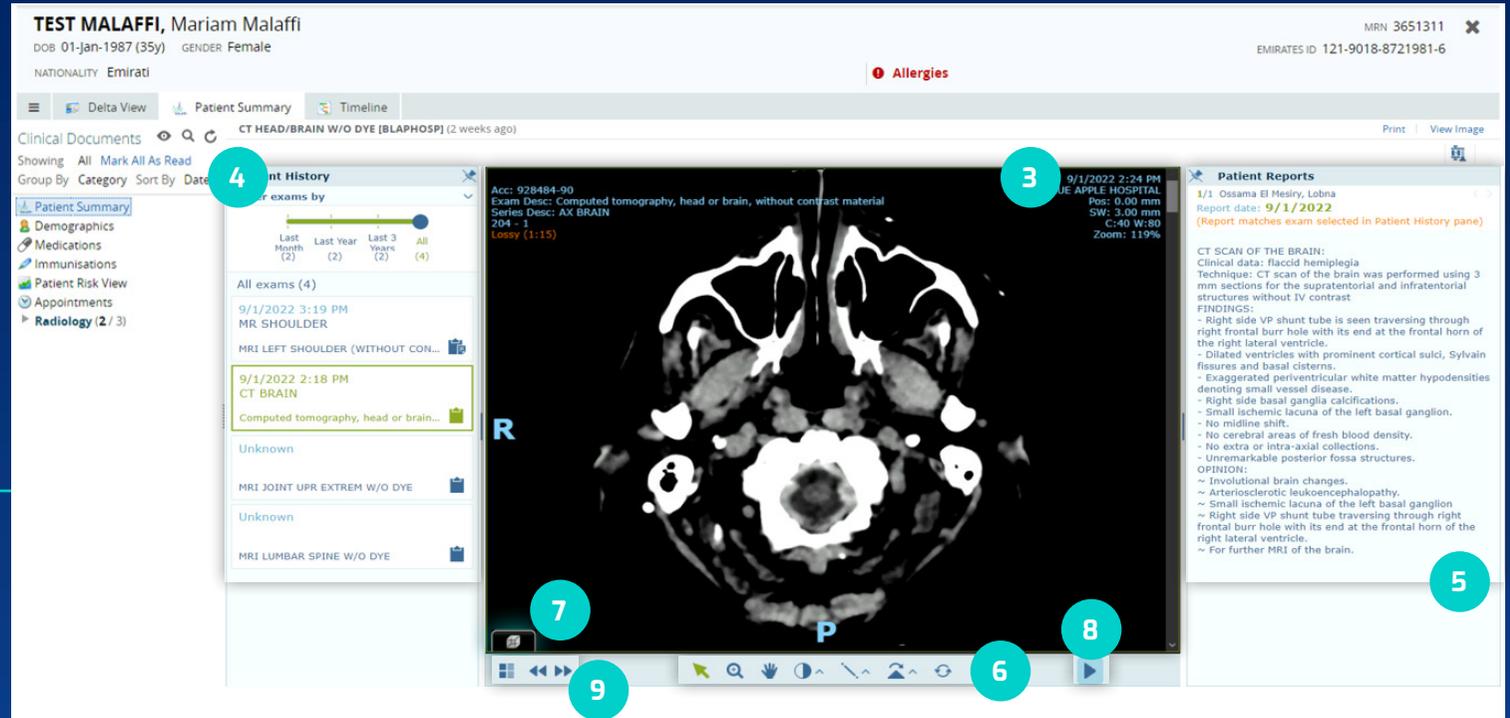

3 Image viewer in iFrame

The image will appear embedded in the Malaffi Provider Portal frame. The main image will display the patient demographic and exam details.

4 Patient History pane

To load a prior exam, open the left pane and click the required exam.

To display a prior exam side-by-side with the current exam, click the + sign that appears when hovering over the top-right corner of the exam box.

If there are multiple historical exams for the patient, these can be filtered by time frame.

To maintain the Patient History pane view, select the Pin icon. To hide, unselect the Pin icon.

5 Patient Reports pane

The right pane shows the report of the displayed exam. To maintain the Patient Reports pane, select the Pin icon. To hide, unselect the Pin icon.

6 Manipulation tools

At the bottom of the Patient Exam screen, there is a toolbar to perform certain operations. The icons that appear here will vary depending on the type of image and the mode of view currently displayed. Hover over the icon to view function or refer to the training guide for more details.

7 Viewing images in 3D

Click the Rendition Options icon in the 3D Tools menu to display different rendition types (MPR, coronal, sagittal, and axial) or pre-defined volume rendering (VoR) protocols.

Note:
3D Rendition tools are only available for CT and MRI Study types.

8 Cine tools

Cine tools are used to control the image motion.

9 Series Chooser

To see a list of all the series available of the displayed exam, the user can select the Series chooser . Click the series that is required to be displayed or drag it to the display area.

Use the arrow icons to navigate to the next or previous series.